

# Agilium® Softfit

Less Pain. More Life.

Ottobock's Agilium portfolio is designed to support patients with knee osteoarthritis (OA) for everyday use, recreational activities, and sports. The Agilium Softfit (50K90) is a single upright knee orthosis with a textile base body for patients with unicompartmental knee OA.

## Features & Benefits

- Proven three-point pressure system causes direct knee offloading and reduces pain
- Adjustable bars for quick and straightforward fitting to treat both varus and valgus types
- Comfortable, slim design and patented vector grip material ensures firm fit without slipping
- Easy to don and doff by the patient themselves


## Indications

Medial or lateral unicompartmental knee osteoarthritis.



## Adjustments

The corrective force can be set with an Allen wrench above the joint to easily adjust the varus/valgus angle of the orthosis.

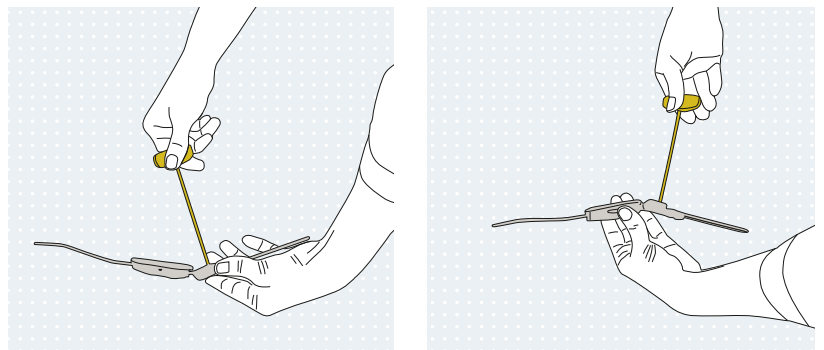

## Ordering Information

| Article number | Side | Size | Circumference |        |
|----------------|------|------|---------------|--------|
|                |      |      | A (cm)        | B (cm) |
| 50K90          | L/R  | S    | 35-38         | 44-48  |
|                | L/R  | M    | 38-41         | 48-52  |
|                | L/R  | L    | 41-44         | 52-56  |
|                | L/R  | XL   | 44-48         | 56-61  |
|                | L/R  | XXL  | 48-51         | 61-67  |

Example format: 50K90=L-S or 50K90=R-M