

VAKC84 | KOOLFLEX™ PULL-UP HINGED KNEE ORTHOSIS

Now available directly from Frontline Medical Services

Previously on VA Denver Logistics Center (DLC) Contract

No: VA791-12-D-0008

FEATURES

- 1/8" Koolflex™ material - breathable, latex-free, hypo-allergenic
- Slip-on design with pull tabs, hook and loop closure
- 2" hook and loop circumferential straps
- Covered aluminum, polycentric hinges
- Reinforced open patella

INDICATIONS

- Mild Anterior & Posterior Cruciate Ligament Instability (ACL, PCL)
- Mild Medial & Lateral Collateral Ligament Instability (MCL, LCL)
- Menisci and Patellar Retinaculum Injuries and Instabilities

SIZING GUIDE

Measure knee circumference with leg extended.

SIZE	KNEE CENTER	6" ABOVE KNEE
XSMALL	12" - 13"	16" - 19"
SMALL	13" - 14"	17" - 20"
MEDIUM	14" - 15"	18" - 21"
LARGE	15" - 16"	19" - 22"
XLARGE	16" - 18"	20" - 24"
2XLARGE	18" - 20"	22" - 26"
3XLARGE	20" - 22"	24" - 28"
4XLARGE	22" - 24"	26" - 30"
CUSTOM	24" and Up	Provide Measurements

CONTRAINDICATIONS

DO NOT wear if you have circulatory problems or are susceptible to dermatitis. If rash, blisters, or skin irritations occur, consult your physician and discontinue use of the support.

(Rev 10/18/2018)



(L-1810)

Made in the USA by New Options Sports

WASHING INSTRUCTIONS: Remove hinges from pockets. Close all hook and loop closures. Hand wash in cold water with mild detergent. NO chlorine bleach. Lay flat and air dry. Do NOT put in dryer.

WARRANTY: Defects in manufacturing or materials used will be covered for a period of (90) days from date of purchase. Items covered by warranty will either be replaced or repaired by the distributor. Supports must be returned CLEAN & DRY. Neither the manufacturer nor distributor will be liable for any damages or injuries arising from defective materials, workmanship, or problems covered under the Indications and Contraindications.

ORDERS: fax 720-221-0677 | orders@frontlinemedicalmil.com | phone 866-681-0444



Frontline Medical Services
Toll Free: 866-681-0444
www.frontlinemedicalmil.com



Serve proud. Live well.™



VAKC84 FITTING INSTRUCTIONS

1. Center support in front of leg and pull straight up so that knee cap is centered in opening.

2. Wrap straps around leg

3. Thread straps through buckles and secure Velcro closure

4. Support is now properly fit

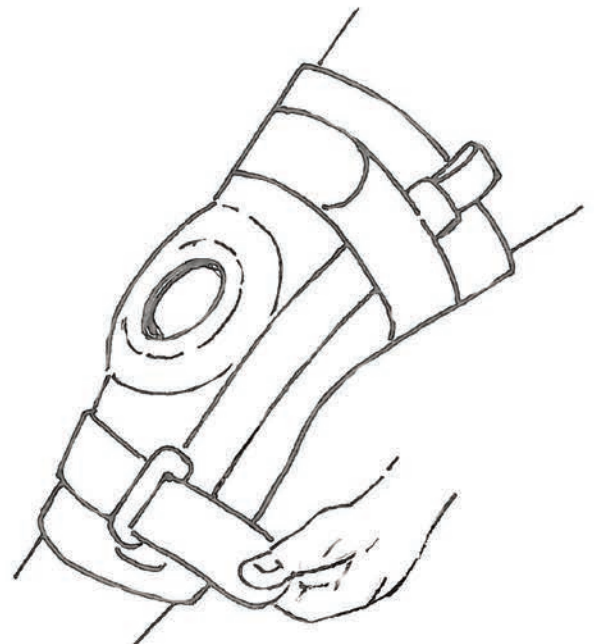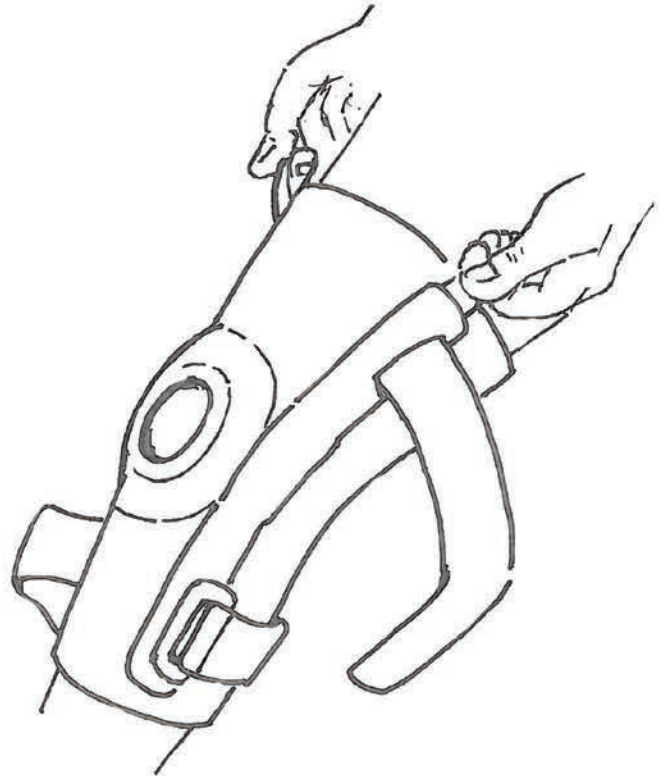

New Options Sports



Made in the USA by New Options Sports

Frontline
MEDICAL SERVICES

Serve proud. Live well.™